




Optokinetic

● Button used during the exercise

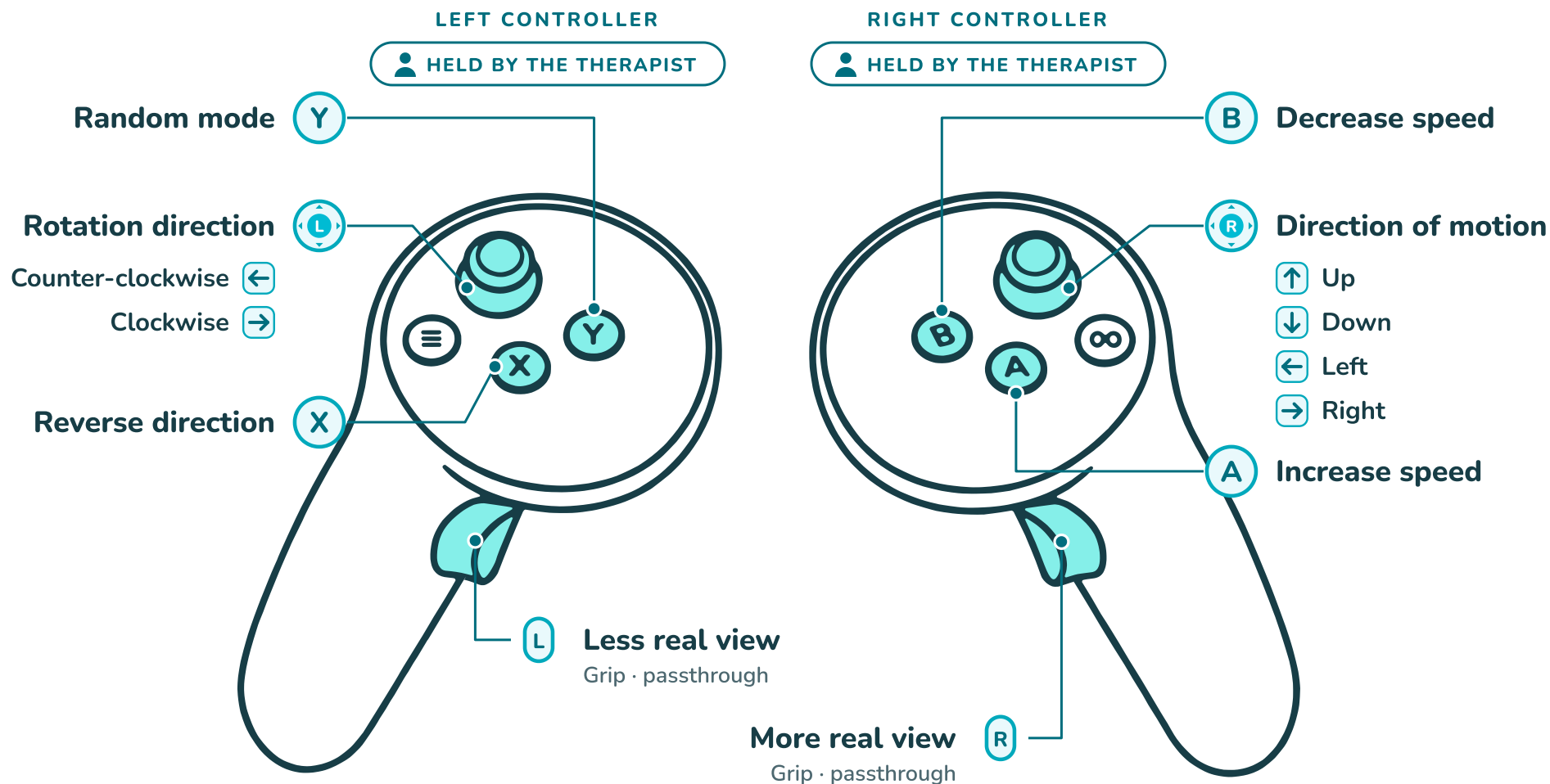

IN EVERY EXERCISE

Right trigger

Short press: pause


Hold 2 s: end the exercise





Translation

● Button used during the exercise

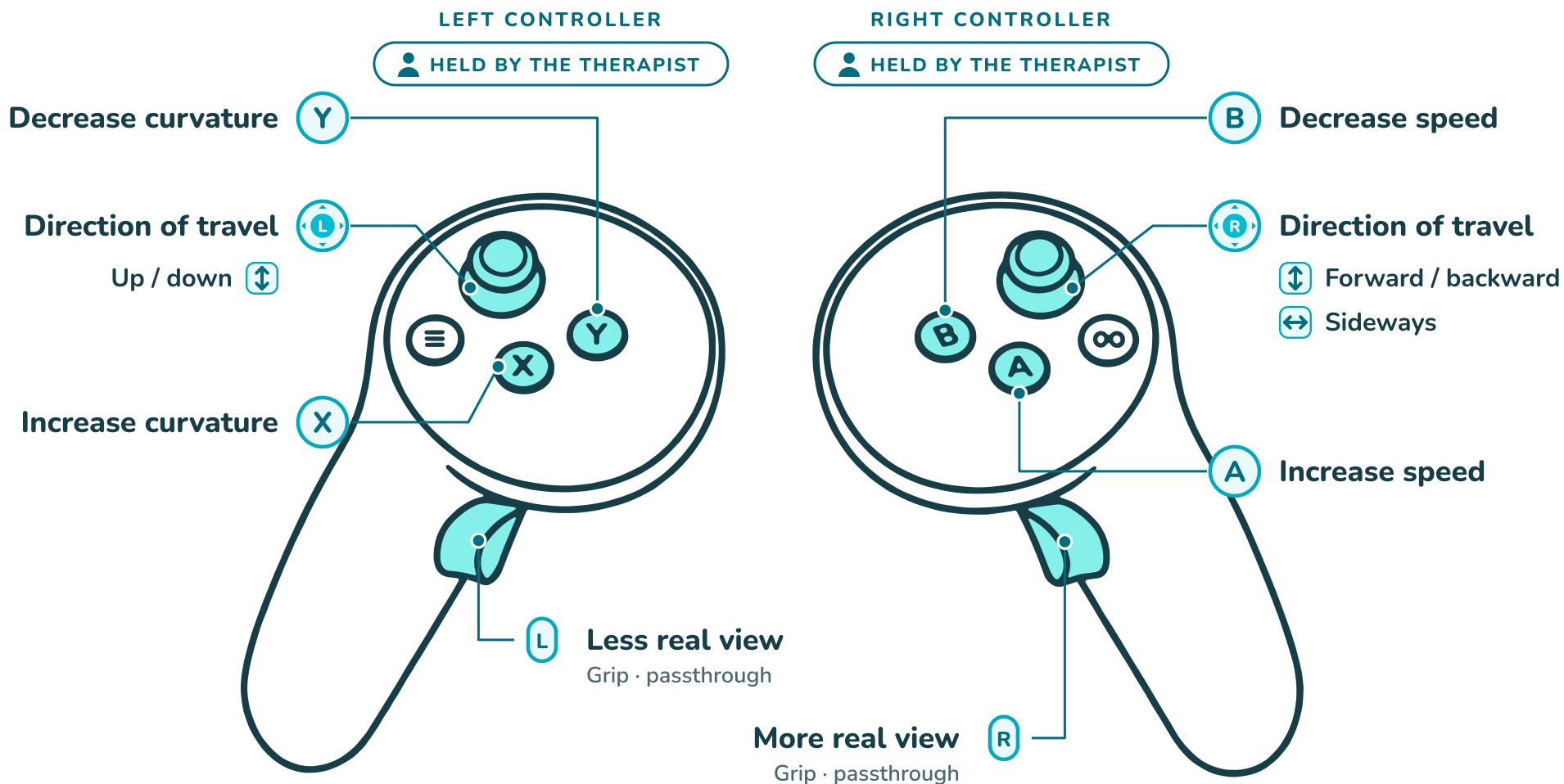

IN EVERY EXERCISE

Right trigger
Short press: pause
Hold 2 s: end the exercise





DVA — gaze stabilisation and shifting

-  Button used during the exercise
-  Controller to hand to the patient


LEFT CONTROLLER

 HAND TO THE PATIENT

Point to the gap
of the C

Up / down 

Left / right 

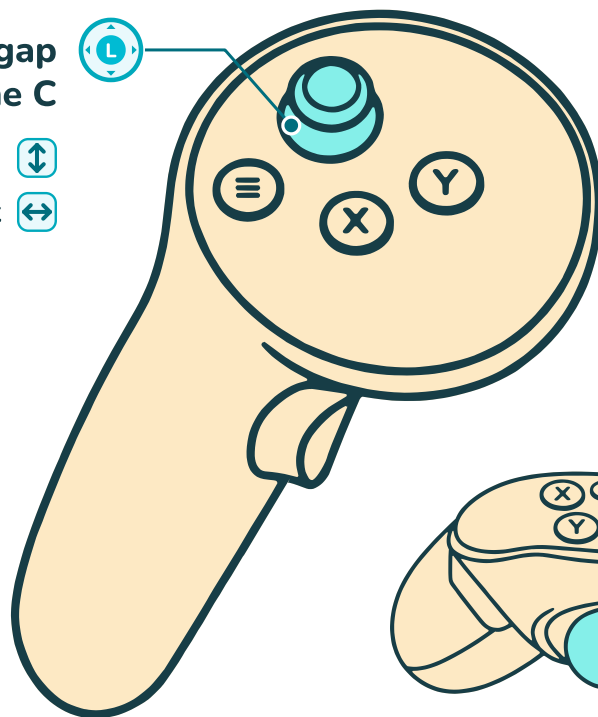


RIGHT CONTROLLER

 HELD BY THE THERAPIST

 B Decrease speed

 A Increase speed

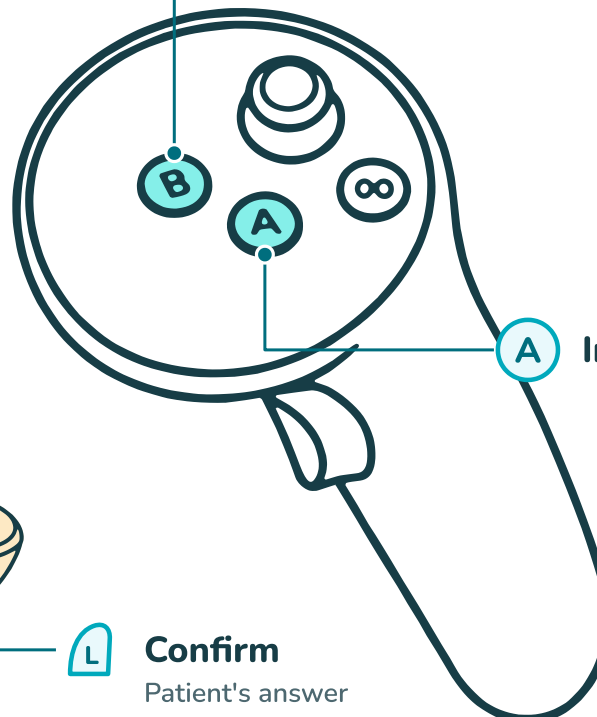

 L Confirm
Patient's answer

IN EVERY EXERCISE

Right trigger

Short press: pause

Hold 2 s: end the exercise




 ∞ **Meta button**
Hold 2 s: recenter the view



Target


● Button used during the exercise

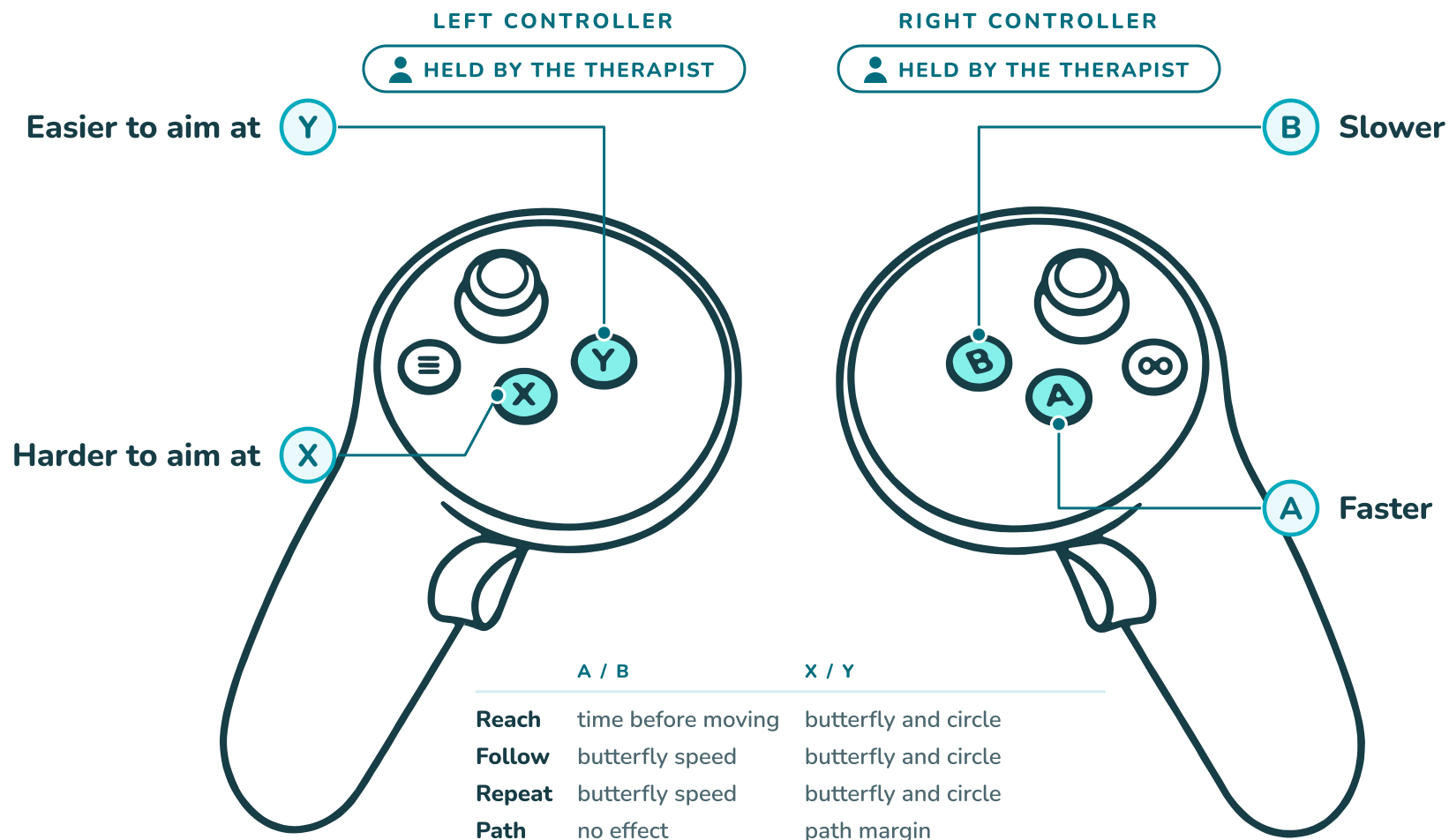

IN EVERY EXERCISE

Right trigger
Short press: pause
Hold 2 s: end the exercise

R

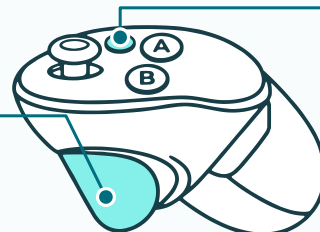

Meta button

Hold 2 s: recenter the view



Visual vertical (SVV)



-  Button used during the exercise
-  Controller to hand to the patient

LEFT CONTROLLER

 HAND TO THE PATIENT

RIGHT CONTROLLER

 HELD BY THE THERAPIST

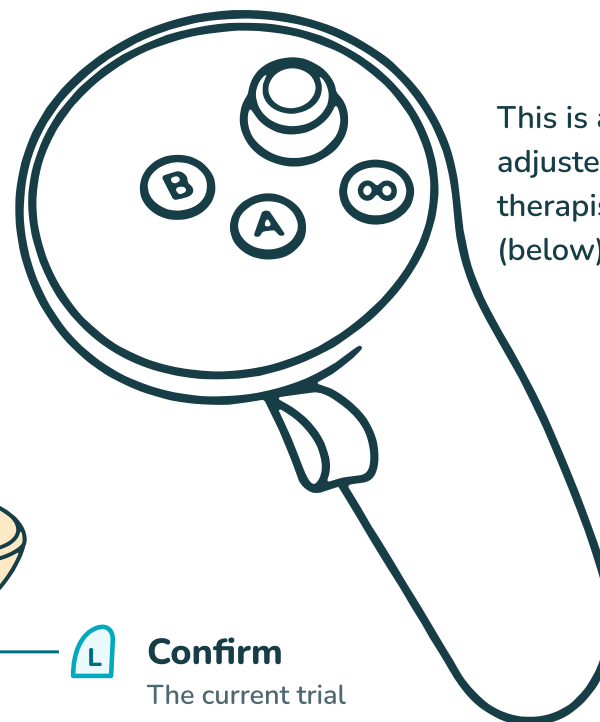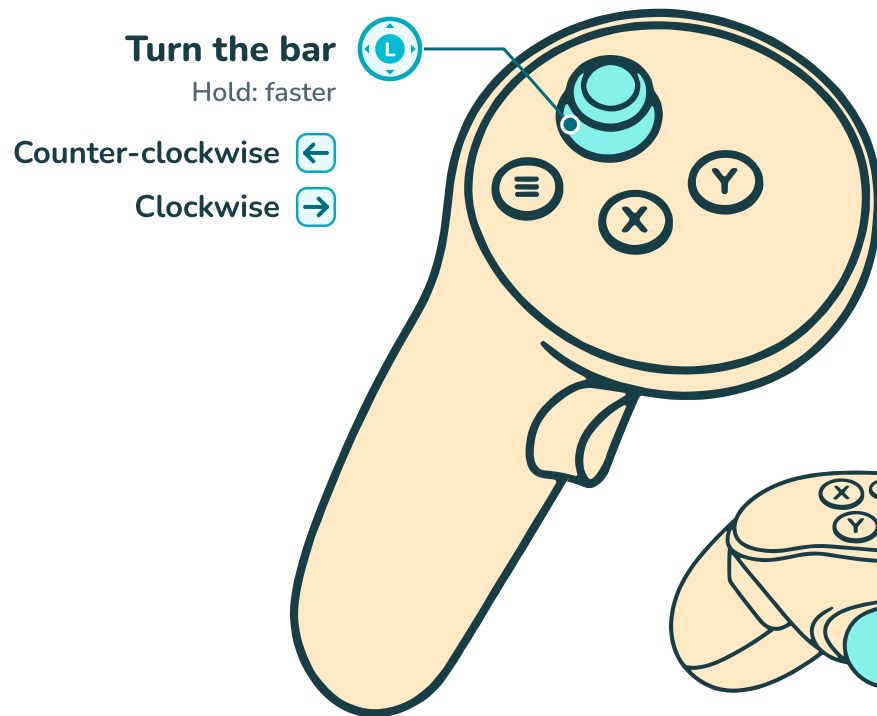

This is a measurement: nothing is adjusted during the series. The therapist only uses the trigger (below).

IN EVERY EXERCISE

Right trigger
Hold 2 s: end the exercise

